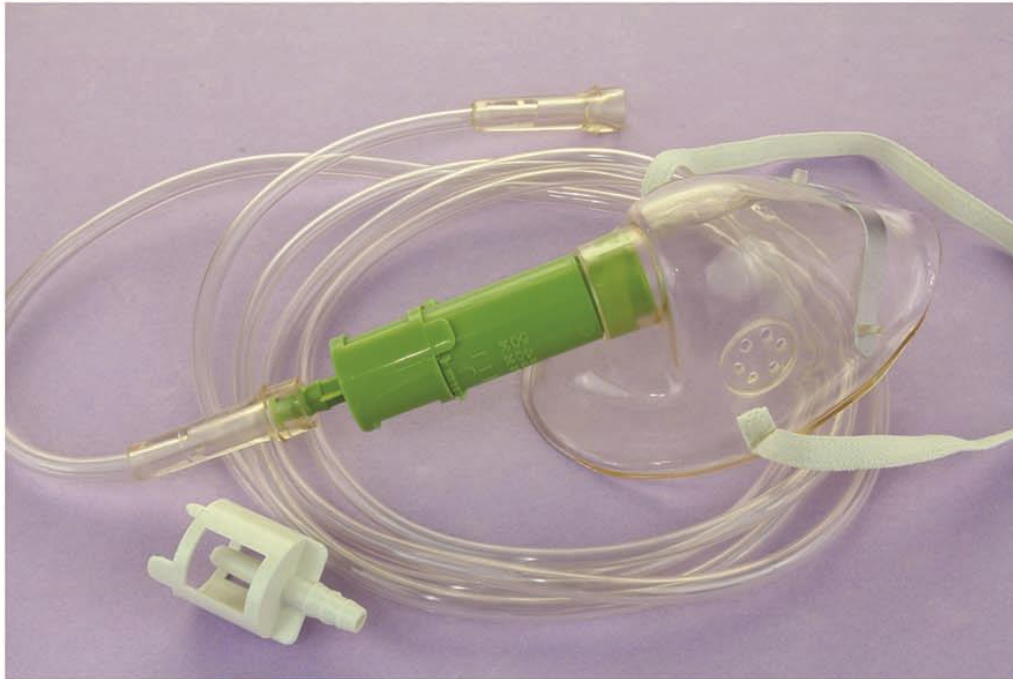


**Ventori Mask
(Oxygen Accurate Flow)**

Code(s) : 1/107

Purpose

Single unit Venturi system that provides five accurate concentrations of oxygen (from 25 to 50 %) " through 2 individual, color coded diluter jets.

Main Features

- Soft, transparent, aerosol mask with anatomical form. Ideal for long term use.
- Wide, adjustable, elastic strap with sealed ends to preclude disassembly of strap.
- Integrated nose clip for proper positioning of mask on the mouth and nasal area.
- Two color- coded diluter jets.
- - Olive-green for 24,26,28 and 30% concentration.
- - White for 35,40 and 50% concentration.
- Easy change of diluter jets in no time.
- Flexible, transparent, anti-c lick tubing of two meter length.
- Both tube ends are fitted with soft funnel-shaped connector for easy connection to oxygen

Sterilization : By EtO gas.